

Retubing Guide


Ear mould

Your Audiologist will have shown you how to retube your ear mould.


This instruction leaflet will support you when you change your own tube.

You will need to change the tubing when it begins to harden and discolour.


The following pictures show you how to change the tubing.

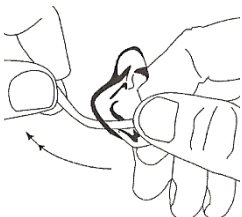


- 1.** Separate the hearing aid and earmould.
Pull the old tubing from the earmould.
Do not throw it away until you have finished.

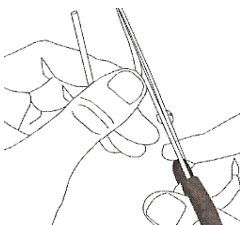


- 2.** Take a new piece of tubing and cut one end of the tube at a slant to make it into a long point.

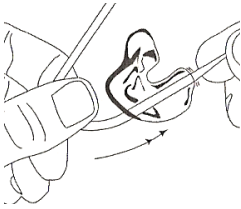


- 3.** Thread the tubing through the earmould. Start from the outside of the mould to the inside. If there are 2 holes thread it into the top hole (which you took the old tube out of).

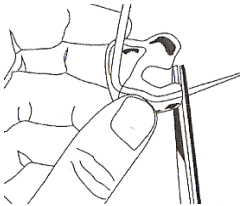

- 4.** Pull the tubing until the other end stands 'upright'.

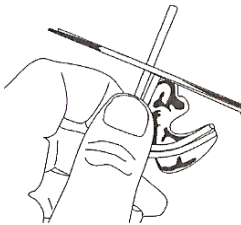

- 5.** Cut the tapered edge of the tubing as close to the earmould as possible. Take the old piece of tubing and line it up with the new tube in the earmould. Cut the new tube to the same length as the old piece.

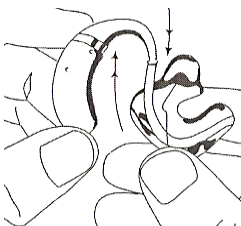

- 6.** Push the end of the new tubing onto the earhook of your hearing aid.

Please visit the Department Website below for a helpful video.

Contact details

DPOW T: 03033 304645

SGH T: 03033 302436

nlg-tr.audiology@nhs.net

www.nlg.nhs.uk/services/audiology

IFP-1238

Date of issue: Nov 20

Date of Review: Nov 23

Author: Ross Keith